



# BeautiCem SA

Fluoride Releasing  
Self-Adhesive Resin Cement



SHOFU



Visit [www.shofu.com](http://www.shofu.com) | Call 800.827.4638

# CASE STUDY



One



Two



Three

Photos courtesy of Shofu, Inc.


BeautiCem SA is a self-etch, self-adhesive, dual-cure resin cement incorporating Shofu's proprietary "Giomers technology". High bond values across all substrates provide piece-of-mind for all of your cementation needs. A quick two-second flash-cure provides easy marginal clean-up. The low 12µm film-thickness allows extra space in tight-fit situations and nearly undetectable margins. The addition of Giomer technology offers well-documented bioactive properties serving to remineralize and protect tooth structure.

## BeautiCem SA

- Self-etching, self-adhesive resin cement
- Unique Giomer technology with the ability to release and recharge fluoride
- No primer required for metal, composite, alumina and zirconia\*
- Auto-mixing syringe eliminates hand mixing
- Dual cure for easy cleanup
- High Radiopacity (1.35 al:mm)
- Low film thickness (<12µm)

Visit [www.shofu.com](http://www.shofu.com) | Call 800.827.4638

Working time (23°/min)	Setting time (37°/min)	Compressive strength MPa	Flexural strength / MPa	
			Light Cure	Chemical Cure
2'50"	3'10"	320	134.4 (10.2)	116.7 (8.0)

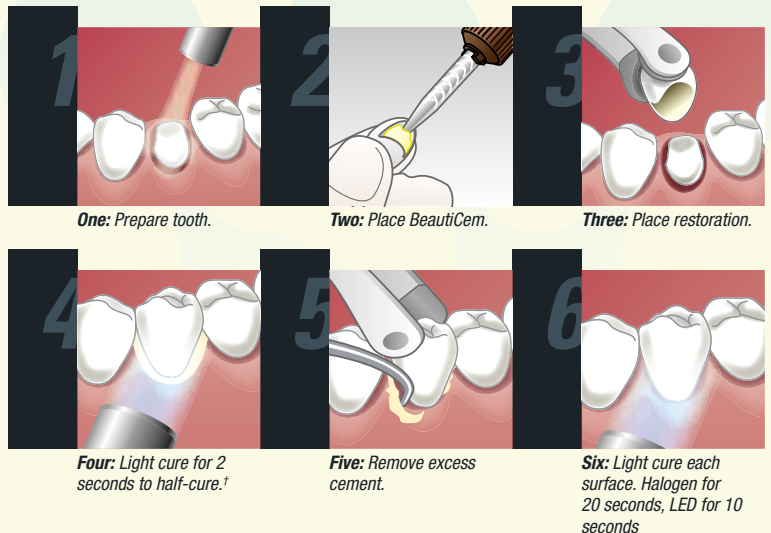

\*Enhanced bonding to all ceramic with Shofu Porcelain Primer.

†For areas where light does not reach, let restoration sit for 3-4 minutes to half-cure chemically.



### Order Information

- Part Number 3213
- Clear
- Part Number 3214
- Ivory
- 5mL paste
- 10 mixer tips
- Step card



Shofu Dental Corporation  
1225 Stone Drive  
San Marcos, CA  
92078-4059, USA

SN33-1018

# Fluoride Releasing Self-Adhesive Resin Cement