

Super-Snap **X-TREME**TM FINISHING & POLISHING & Super-Snap[®]



VISIT www.shofu.com | CALL 800.827.4638

SUPERIOR FINISHING

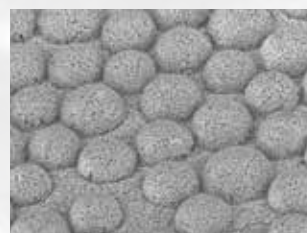


Super-Snap X-TREME FINISHING & POLISHING

Advanced 3D Technology

Super-Snap X-Treme, the latest generation in X-Treme polishing, is designed to x-ceed the aesthetic x-pectations of patients today. Super-Snap X-Treme aluminum oxide disks offer x-tra thickness for greater flexibility and durability, and a new 3D semispherical grit to provide X-Treme gloss results. Super-Snap X-Treme—taking polishing to the next level!

- An enhancement to the original green and red Super-Snap disks
- Added thickness for improved durability and tactile feel while maintaining flexibility interproximally, without warping
- 3D X-Tra coating combined with semi-spherically shaped grits of red disks prevent clogging, secondary scratches and material's denaturation
- Achieve lustrous results in as little as 15 seconds, with a simplified 2-step X-treme procedure⁺
- 2 sizes: standard 12 mm and mini 8mm, double-sided, disposable aluminum-oxide disks
- For a complete 4-step system use with Super-Snap coarse (black) and medium (violet) disks



400x 200µm
3D semispherical coating



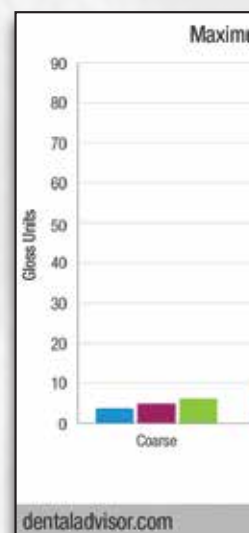
Red Disk
Polished tooth surface x400



Super-Snap X-Treme Disk
Polishing
Photo courtesy of Sam Halabo, DDS



Super-Snap X-Treme Disk
Super-polishing
Photo courtesy of Sam Halabo, DDS




*Test results on file. All trademarks are the property of their respective owners.
Super-Snap X-Treme pro
Sof-Lex Spirals



Super-Snap X-Treme
Technique Video

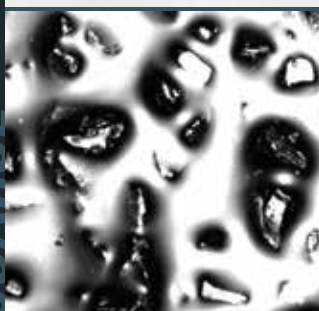
COARSE



Super-Snap Black Disk
Coarse


Laser microscope disk surfaces (x400)

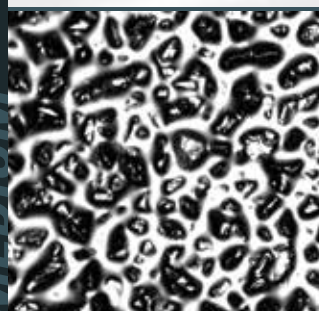MEDIUM




Super-Snap Violet Disk
Medium

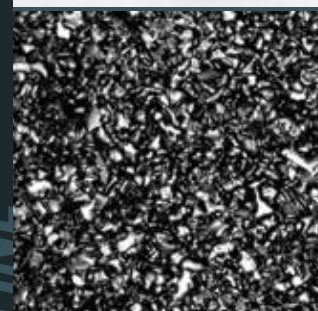FINE




Super-Snap Green Disk
Fine

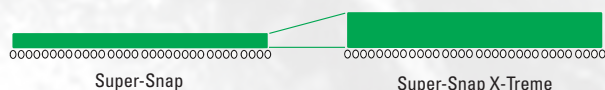& POLISHING DISKS




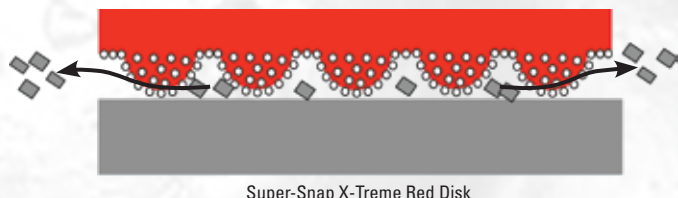Super-Snap X-Treme polishing disks are thicker than Super-Snap allowing for greater force to be applied during polishing. The added thickness creates a greater tactile feel and produces higher gloss results.



Super-Snap X-Treme's red 3D aluminum oxide grits allow space for ground debris to be discharged during polishing. This helps prevent secondary scratches, clogging, and reduce generated heat without any denaturing of the material.

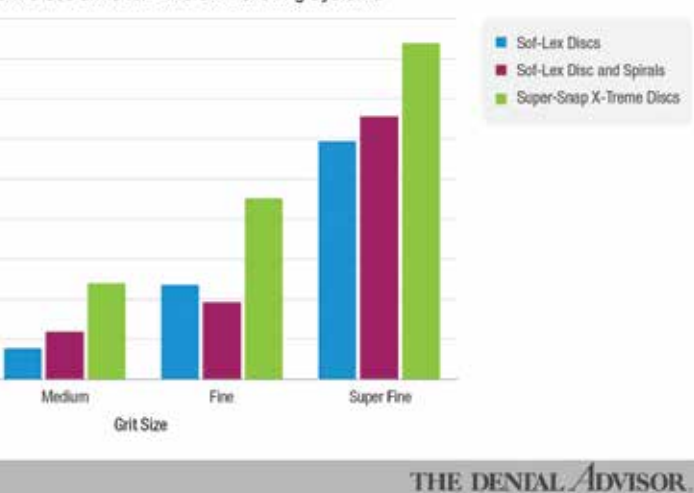Gloss Units for Different Polishing Systems



Results belong to their manufacturers.

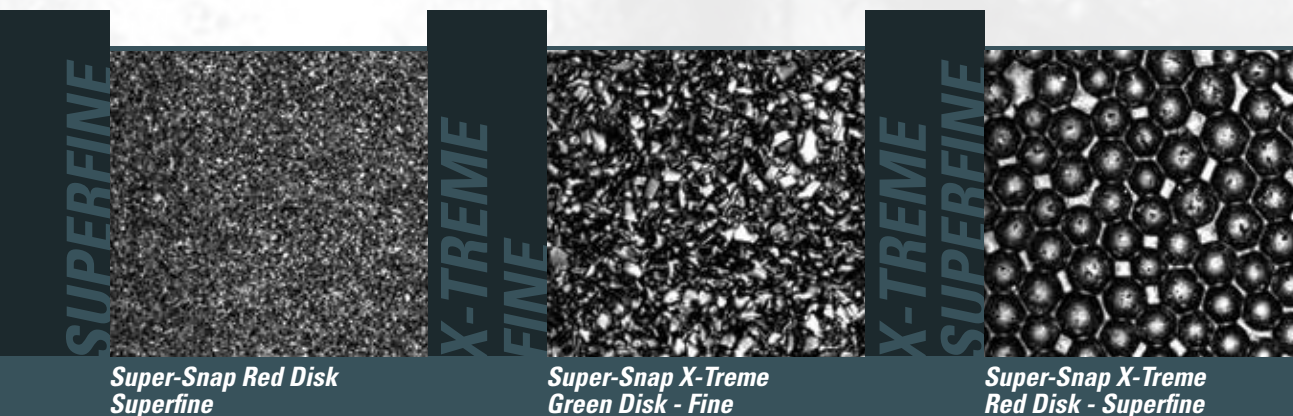
Super-Snap X-Treme produced 41% greater gloss than Sof-Lex Discs and 29% greater gloss than

Super-Snap











"Top Composite Finishing and Polishing System" Awarded by THE DENTAL ADVISOR for 6 consecutive years

Super-Snap, the 4-step disk system, is designed for easy, fast and safe contouring, finishing and polishing of all microfilled and nano-hybrid composites. The grit consists of silicon carbide or aluminum oxide, without a metal center to risk gouging. Super-Snap disks are ensured to deliver superior polish and high gloss results. Super-Snap Polystrips offer the same results as the Super-Snap disks but are double-ended strips, with a different grit at each end for interproximal areas.

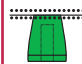

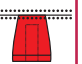

- Easy to use 4-step contouring, finishing and polishing system
- Color coded by grit size, course (black), medium (violet), fine (green) and superfine (red)
- Ultra-thin, high flexibility is ideal for interproximal use
- No metal center eliminates risk of gouging
- Elastic silicone shank mount for easy placement on the mandrel
- Disposable double sided disks in 2 sizes – 4 grits (silicon carbide and aluminum carbide)
- Available as a kit or single packs



Super-Snap
Technique Video

SHAPE CATEGORY		Contouring		Finishing				Polishing		Super-Polishing	
Super-Snap Shape											
Color/Grit		Black (coarse)		Violet (medium)				Green (fine)		Red (super fine)	
Refills (50 pk)	PN	L506	L507	L508	L528	L509	L522	L501	L519	L502	L521



SHAPE CATEGORY		Polishing		Super-Polishing	
Super-Snap X-Treme					
Shape					
Color/Grit		Green (fine)		Red (super fine)	
Refills (50 pk)	PN	L511	L531	L512	L532

VISIT www.shofu.com | CALL 800.827.4638

Order Information

Super-Snap X-Treme Kits

- Part Number 0506 – 50 green standard X-Treme disks, 50 red standard X-Treme disks
- Part Number 0507 – 50 green mini X-Treme disks, 50 red mini X-Treme disks

Super-Snap X-Treme Ultra-Gloss Performance Kits

- Part Number 0500XS – 100 Super-Snap standard disks, 100 Super-Snap X-Treme standard disks, 25 Super-Snap Buff disks, 4 CA metal mandrels, DirectDia polishing paste, **BONUS:** F00 A2 Beautifil Flow Plus and compartment case
- Part Number 0500XM – 100 Super-Snap mini disks, 100 Super-Snap X-Treme mini disks, 25 Super-Snap Buff disks, 4 CA metal mandrels, DirectDia polishing paste, **BONUS:** F03 A2 Beautifil Flow Plus and compartment case

Super-Snap Kits

- Part Number 0500 – Rainbow Technique Kit 100 standard disks, 80 mini disks, 40 polystrips, 2 Dura-White Stones (CN1 and FL2),

CompoSite fine midi-point, 4 CA metal mandrels

- Part Number 0505 – Mini Kit 32 standard disks, 16 mini disks, CompoSite fine midi-point, CA metal mandrel

Super-Snap Singles – Individual Patient Use

- Part Number 0501 – Singles 2 disks of coarse, medium, fine and super fine grits, Super-Snap plastic CA mandrel
- Part Number 0502 – Standard 1 disk of coarse, medium, fine, super fine grits, CA mandrel
- Part Number 0503 – Mini 1 disk of coarse, medium, fine, super fine grits, CA mandrel
- See chart above for individual shapes

Super-Snap Buff and Super-Snap SuperBuff Kits

- Part Number L523 – Super-Snap Buff Disk 25 disks
- Part Number L524 – Super-Snap Buff Mini Disk 25 disks
- Part Number 0535 – SuperBuff Set 25 disks and 2 CA mandrels
- Part Number 0536 – SuperBuff Disk 50 disks (refills)

Super-Snap Polystrips

- Part Number L525 – Polystrips 100 black (coarse) and violet (medium)
- Part Number L526 – Polystrips 100 green (fine) and red (super fine)

Mandrels

- Part Number 0439 – CA Metal
- Part Number 0440 – CA Plastic



Super-Snap Disks

+In a Dental Advisor bench test, all direct restorative materials polished with Super-Snap X-Treme disks demonstrated an above-average sheen within just 15 seconds. The attained gloss values were about 20% higher as compared to a 4-step polishing system and nearly 3-times better than those accomplished with a one-step procedure.

USA

Shofu Dental Corporation
1225 Stone Drive
San Marcos, CA 92078-4059, USA

JAPAN

Shofu Inc.
11 Kamitakamatsu-cho
Fukuine, Higashiyama-ku
Kyoto 605-0983, Japan

SINGAPORE

Shofu Inc. Singapore Branch
10 Science Park Road #03-12
The Alpha, Science Park II
Singapore 117684

U.K.

SHOFU UK Riverside House,
River Lawn Road,
Tonbridge,
Kent, TN9 1EP, UK

CHINA

Shofu Dental Supplies (Shanghai) Co. Ltd.
No.645 Jiye Road
Seshan Industry Park
Songjiang 201602, Shanghai, China

GERMANY

EU Representative:
Shofu Dental GmbH
Am Brüll 17
40878 Ratingen, Germany